

# Water Purifier - Auto Pure Type 1 Water

For General Chemistry, Buffer Solutions, Academia



**Benchtop**                      **WA301B: 120V / WA311B: 220V**  
**Remote Dispense**       **WA301R: 120V / WA311R: 220V**

Water quality	Type 1, 18.2 MΩ	Feed requirement	< 20 uS	Flow	0.67 gpm 2.5 lpm
---------------	-----------------	------------------	---------	------	---------------------

Water is purified using a stage purification process which includes high-purity ion exchange resins to remove dissolved minerals and internal recirculation to maintain purity. As water exits the system, a final filter removes particulates and bacteria to reach Type I water specifications.



### ■ WA Series Configuration

- Conditioning cartridge
- Polishing cartridge
- 0.2 micron filter

### ■ Features

- 2.5 LPM\* of 18 MΩ Type 1 reagent water on demand
- Intuitive, user-friendly touch screen interface
- Real time quality monitoring with leak detection
- **Compact design with remote dispense option**
- **Quick change, no tools cartridge design**
- Large capacity, low cost consumable cartridges
- Premium quality ion exchange resin
- **Low ownership cost**

\* Reported flow rate is typical but can vary depending on supply pressure and system options

### ■ Specifications

Model	WA301B/WA311B	WA301R/WA311R
Type	Benchtop	Remote dispense
Resistivity*	18.2 MΩ-cm	
Bacteria*	< 1 cfu/ml	
Particulates*	< 0.2 μm filtration	
TOC*	< 15 ppb	
Temperature	100°F / 30°C	
Pressure	90 PSIG Max. / 20 PSIG Min.	
Flow rate (typical)	0.67 Gallons Per Minute / 2.5 Liters Per Minute < 2.0 LPM with Ultrafilter & Endotoxin filter installed	
Inlet connection	3/8" O.D. Tubing / 9.53 mm O.D. Tubing	
Dimensions (H x W x D)	21.6 in. x 14.7 in. x 14.4 in. (54.8 cm x 37 cm x 36.5 cm)	
Weight	42 lbs. dry / 47 lbs. operating (22.7 kg dry / 25 kg operating)	
Power requirements	120V, 60 Hz @ 1.0 amp / 220V, 50 Hz @ 0.5 amp	


\* Performance based upon RO / DI or Distillation Feed (< 20 uS/cm, < 50 ppb of TOC, 0.2 micron filtration)

### System Highlights


Intuitive Touch Screen Controls  
 Programmable Volume Dispense  
 Real-time Quality Monitoring & Alerts

Quick Change Cartridge Design

Large Capacity Cartridges filled with Premium ResinTech Media

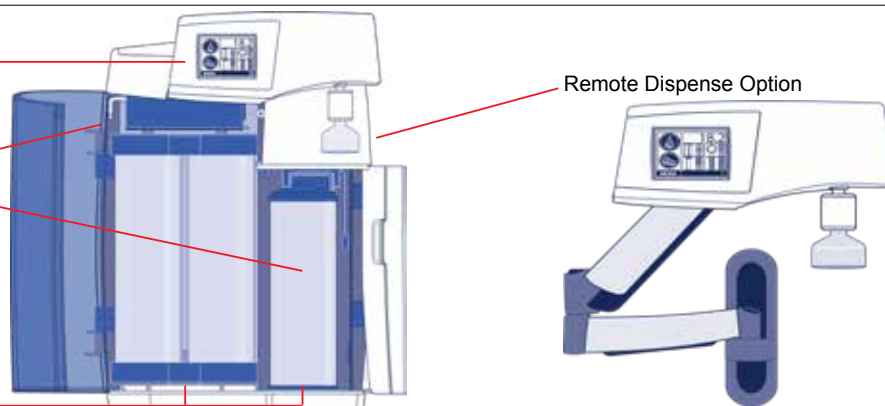

## Flow Diagram

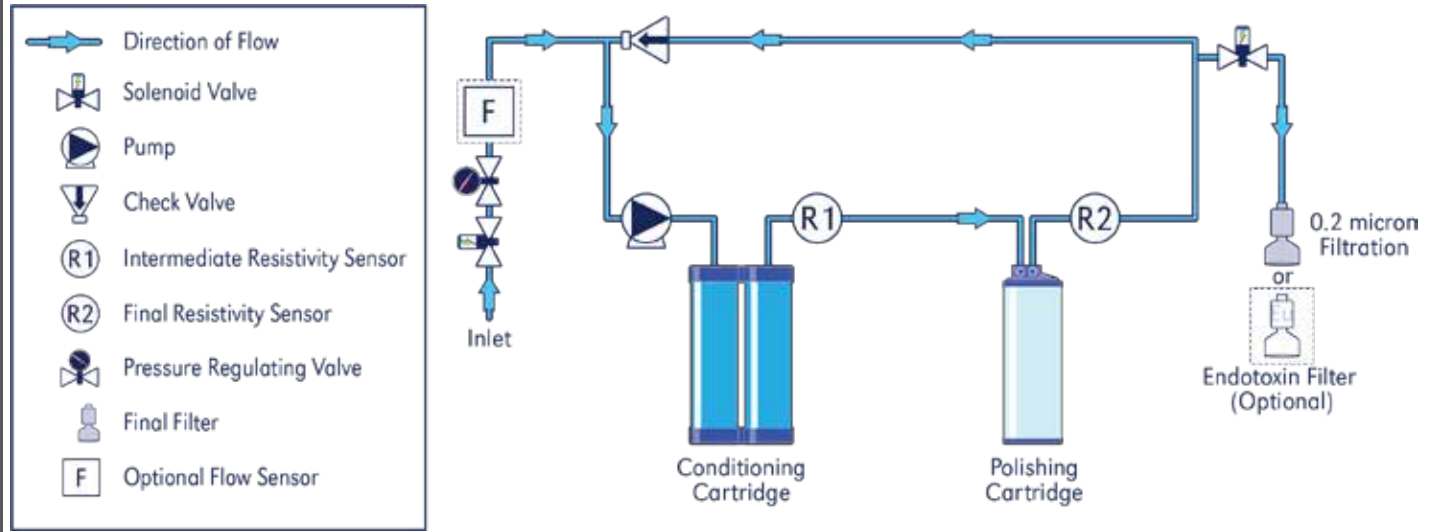

### Additional System Components

Product code	Description
<b>Installed Options and Accessories*</b>	
ARI-PHADF	Direct Feed Port
ARI-PHADG	Recirculating Dispensing Gun
ARI-PHAWB	Wall Mount Bracket
<b>Cartridges and Filters**</b>	
ARI-PX135001	Dual Conditioning Cartridge (RO Feed)
ARI-PX135002	Dual Conditioning Cartridge (SDI / TAP Feed)
ARI-PX115103	Polishing Ultrapure Cartridge
ARI-PX115104	Polishing Ultrapure Low Organics Cartridge
ARI-PF006402	0.2 Micron Capsule Final Filter - 1/4" MNPT x 1/4" hose barb
ARI-PF007101	0.2 Micron Capsule Final Endotoxin Filter - 1/4" MNPT x 1/4" hose barb
<b>Consumables</b>	
ARI-PF007105	Ultrafilter - 3/8" Quick Disconnect Fitting
ARI-HPA016	UV Bulb 254/185nm TOC Destruct
ARI-HPA017	Sanitization Kit (includes cartridges, test strips, and gloves)

\* Must be pre-ordered for factory installation

\*\* Not included in a system and ordered separately

### Best Value in Replacement Consumables

